



Internal Organ Soft and Flexible Retractor

Silicone Disc

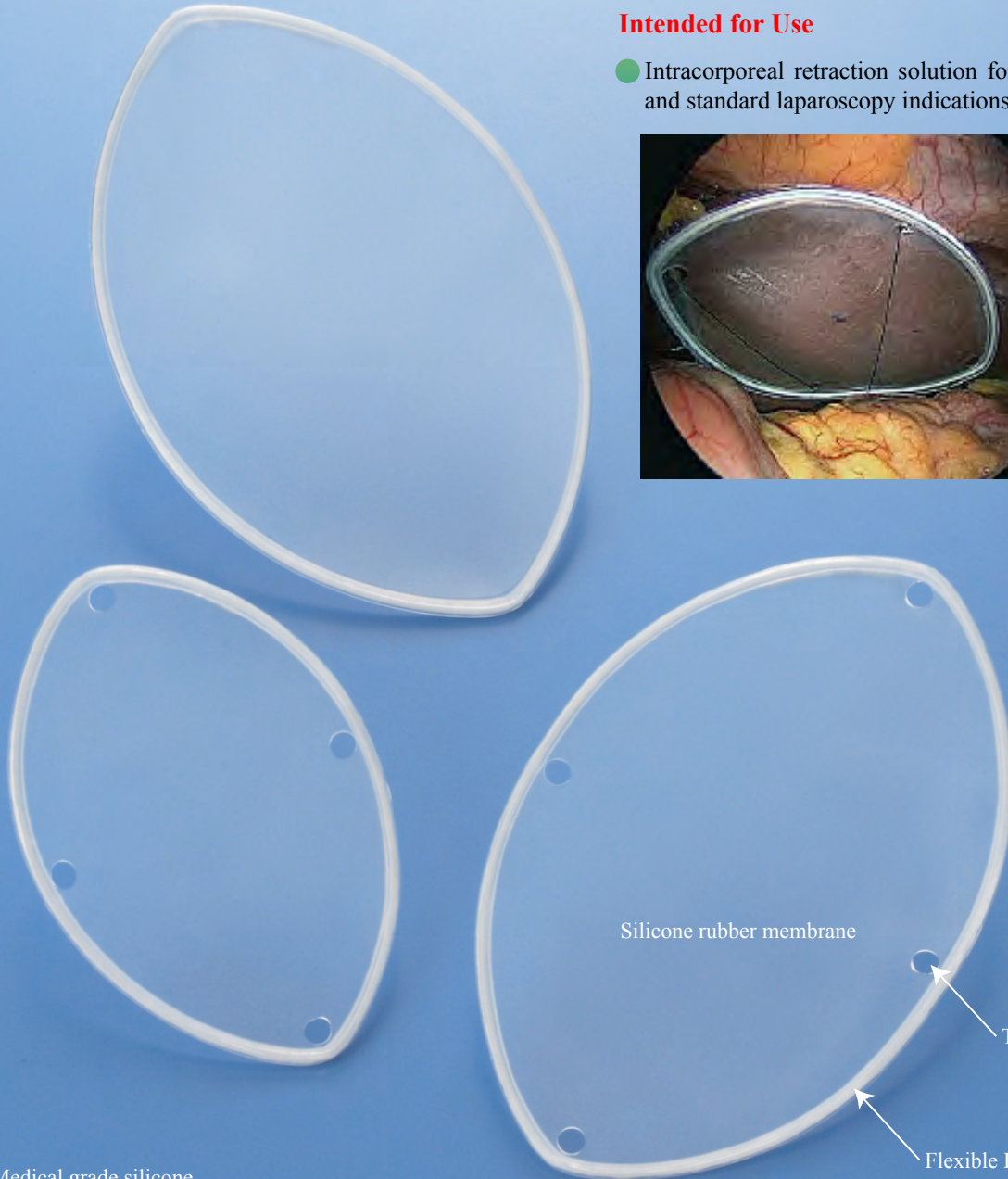

Intended for Use

- Intracorporeal retraction solution for single incision and standard laparoscopy indications.

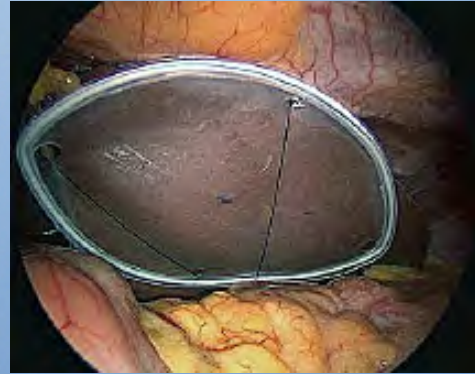

Material: Medical grade silicone

Features

- This consists of a silicone rubber membrane inside a flexible ring and 4 holes for thread lifting.
- The disc can be inserted into the abdomen under the pneumoperitoneum through a 12mm port or small incision(2-3cm) space.
- An organ retraction during laparoscopic upper abdominal surgery would lift the organ out of the operative field in a nontraumatic manner, providing adequate space for visualization and operation for as long as required, and without the need for additional ports or assistance.

Specification	
Device Code	Description
29111400	Elliptical 160x100mm
29111410	Elliptical with 4 holes(Large) 160x100mm
29111420	Elliptical with 4 holes(Small) 120x80mm
*29111450	Heart-shaped with 3 holes 100mmx120mm
Packaging	
2units per box, individual peel pack for single-use with Eto sterilization.	

*)Order made only

Please read the instructions and all accompanying documents carefully prior to use to avoid any failure.

Distributed by:

Manufacturing site;
Hakko co.,Ltd.
1490 Isobe, chikuma-shi, Nagano, Japan
www.hakko-medical.com

